
miniland

baby tooth care

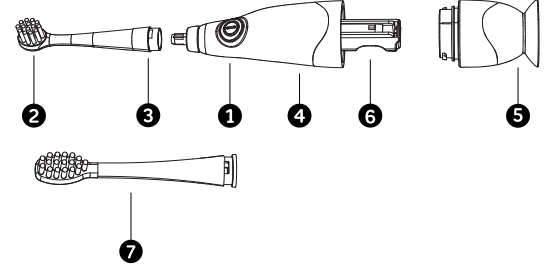ESPAÑOL

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ENGLISH

1. SAFETY PRECAUTIONS

Whenever using any electrical device, it is essential to follow basic safety precautions, including the following:

- READ ALL SAFETY INSTRUCTIONS CAREFULLY BEFORE USING THE PRODUCT.
- THE FEATURES DESCRIBED IN THIS MANUAL ARE SUBJECT TO CHANGE WITHOUT PRIOR NOTICE.**
- Designed for children from their first tooth. Must be used under adult supervision. This is not a toy.
- Children and individuals with reduced abilities should use it under supervision. Assistance is recommended for younger children.
- Do not use if damaged. Inspect and replace the brush head every 3 months or sooner if worn out.
- Do not modify or disassemble the product, as this may compromise safety and void the warranty.
- Water-resistant (IPX7): Can be rinsed but should not be submerged or exposed to corrosive liquids.
- Do not use while the child is walking or running, nor on wounds or inflamed gums without a dentist’s recommendation.
- Keep out of reach of young children when not in use.
- Use only AAA alkaline batteries, ensuring correct polarity. Do not mix new and used batteries or use damaged batteries.
- Remove the batteries if not in use for an extended period and dispose of them properly.
- Do not use in the bathtub or shower if there is a risk of falling into deep water.
- Avoid extreme temperatures, heat, or fire.
- To clean, rinse with water and dry with a towel. Do not use harsh cleaning products.
- If damaged or malfunctioning, discontinue use and contact customer service.

WARNING!

- Keep the device dry, clean, and away from heat.
- If dropped, inspect the brush head before using it.
- Consult your dentist if undergoing oral treatment.
- Store the toothbrush and its accessories out of reach of young children.
- Personal use only; not intended for dental clinics or institutions.

2. DESCRIPTION, COMPONENTS, AND CONTENTS

- Power button
- Small round brush head with a white ring (From the first tooth up to 24 months)
- Light
- Soft-grip handle

- Suction cup
- Battery compartment
- Large elongated brush head (From 24 months to 6 years)

The box also contains:

- 1 AAA battery, 1 instruction manual, and warranty

If any of the above-listed contents are missing, please contact your point of purchase.

3. INSTRUCTIONS FOR USE

3.1. ABOUT THE BRUSH HEADS

The product includes two brush heads designed to suit your child’s age. Each head has a different shape and size, with a distinctive color ring (green or white) for easy identification:

- **Small round brush head** (2): For children from their first tooth up to 24 months. The toothbrush usually comes with this head pre-installed.
- **Large elongated brush head** (7): For children from 24 months to 6 years.

The toothbrush features a water-resistant handle (4), is electrically safe, and is designed for use in the bathroom.

3.1.1. ATTACHING/REPLACING THE BRUSH HEAD

An adult should perform this process by following these steps:

To attach the brush head:

- Hold the brush head (2 or 7) and the handle (4) without touching the bristles.
 - Turn the handle (4) to the left so the brush head (2 or 7) rotates 45° to the right.
 - Turn the handle (4) to the right until the brush head (2 or 7) is aligned.
- To replace the brush head:**
- Hold the brush head (2 or 7) and turn the handle (4) to the left.
 - Remove the brush head (2 or 7) and insert the new one.
 - Turn the handle (4) to the right until it is properly aligned.

NOTE 1: Ensure the brush’s notch aligns with the brush head’s notch for a secure fit.

NOTE 2: Replace the brush head every 3 months or sooner if worn.

NOTE 3: Miniland reference: 89679 “Replacement brush heads”.

3.2. OPERATING MODES

3.2.1. LIGHT-ONLY MODE

The toothbrush features a light mode to check brushing quality.

- Press the power button (1) once to activate it.
- The light (3) will stay on for 30 seconds and then turn off automatically.

3.2.2. TIMER BRUSHING MODE

The Miniland toothbrush includes a 2-minute timer to ensure proper brushing.

- Moisten the bristles, apply a small amount of child-friendly toothpaste, and place the brush head on the child’s teeth.
- If the toothbrush is off, press the button (1) twice. If light mode is active, press it once.
- Brush the lower teeth first (1 min), then the upper teeth (1 min), covering all surfaces.
- The toothbrush will automatically turn off after 2 minutes, ensuring the recommended brushing time. To stop it earlier, press the power button (1).

NOTE: Avoid pressing too hard. If the bristles wear out before 3 months, excessive pressure may be applied. Light bleeding at the beginning is normal and should subside within a few days. If it persists for more than two weeks, consult a dentist.

4. MAINTENANCE

Usage and Cleaning

After each use:

- Rinse the brush head under running water with the handle turned on for a few seconds.
 - Turn off the toothbrush and remove the brush head, following the steps in section 3.1.1 of this manual.
 - Wash the brush head and handle separately under running water.
 - Dry them well and reassemble both parts.
- The toothbrush features a suction cup (5) at the base for convenient bathroom placement. Stick it to a smooth, dry, and non-porous surface for better adhesion.
- Battery Replacement**
- The battery compartment (6) is located at the base of the toothbrush.
- Firmly hold the toothbrush and turn the bottom part to the left.
 - Insert one AAA alkaline battery (1.5V), ensuring correct polarity.
 - Close the lid to prevent water from entering.

5. TROUBLESHOOTING

- The device does not turn on: Check battery placement and replace if necessary.
- The toothbrush turns off after 30 seconds: This is normal when in light-only mode.
- The toothbrush turns off automatically after 2 minutes: This is normal.
- For any other issues, contact Miniland customer service.

6. BATTERY AND PRODUCT DISPOSAL INFORMATION

Before disposing of the electric toothbrush, remove the batteries and dispose of them responsibly. Do not mix products marked with the crossed-out waste bin symbol with household waste. For proper disposal, take them to an authorized collection point or consult your point of purchase.

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