miniLAND
BABY

thermotalk plus

Manual de Instrucciones • Instruction Manual • Manual de instruções
Mode d'emploi • Bedienungsanleitung • Manuale d'istruzioni
Instrukcja obsługi • Návod k obsluze • Руководство по эксплуатации

ÍNDICE

ESPAÑOL

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INDEX

ENGLISH

1. Introduction
2. Contents
3. Safety instructions
4. Product features
5. Instructions for use
6. Maintenance
7. Technical specification
8. Battery & product disposal information

1. INTRODUCTION

Congratulations on your purchase of this clinical IR thermometer with pulse reading function.

Before using this device, please read all instructions in this manual carefully. For the safe and correct use of this product, please read the safety precautions explained in this user manual.

Please store for future reference as it contains important information.

2. CONTENTS

- 1 clinical IR thermometer
- 1 support
- 1 storage and travel bag
- 2 x AAA batteries
- 1 instruction manual and guarantee

If any of the above is missing, please contact your retailer.

3. SAFETY INSTRUCTIONS

3.1. SAFETY PRECAUTIONS

Warning marks and symbols are included for safety and a correct use of this product and prevention of any injury.

See below for the description of warning marks and symbols:

- **WARNING:** Means the possibility of personal injury in case of improper use.

- **NOTICE:** Means the possibility of injuries or damage* in the case of improper use.

*[*This covers any damage to your house, personal belongings and pets]*

• **PROHIBITION:** Means that the detailed items expressed in words or figures within or beside the symbol are totally forbidden. The symbol on the left means forbidden in general.

• **MUST BE OBSERVED:** This means obligatory with the detailed items expressed in words or figures within or beside the mark. The one on the left means that it is compulsory in general.

• **THE MEANING OF THE SYMBOL:** Means Type-B device.

3.2. WARNING

• It is dangerous for patients to perform a self-evaluation and self-treatment based on the results of the reading. Be sure to follow your doctor's instructions.

**A self-evaluation may cause deterioration in the state of the patient's health.*

• Don't touch or blow on the infra-red sensor.

**A dirty infra-red sensor may cause an inaccurate reading.*

• Clean the infra-red sensor with a soft wet cloth.

** Do not clean with toilet paper or a tissue as this may scratch the infrared sensor and cause an inaccurate reading.*

• Keep the machine out of children's reach

** If a child takes his own reading, it may cause damage to his ear. Please consult the doctor at once if a battery is accidentally swallowed.*

• In the case of having performed any outdoor activities, no reading should be taken until the patient's body temperature is the same as the room temperature.

**If the patient has used a water-filled pillow or hot water bottle, this may cause an inaccurate reading.*

• If there is a difference in temperature between the storage area and the place where the reading is to take place, leave the thermometer at room temperature for approximately 30 minutes.

**A failure to do so may cause an inaccurate reading.*

• If the reading exceeds 38°C and / or the patient feels uncomfortable, the user should get in contact with a doctor.

- Stop using the product in the event of any pain.
**It is likely to cause damage to the inside of the patient's ear.*
- It is not recommended to use the product by users suffering ear disease including otitis and tympanitis.
**The affected part is likely to get worse.*
- Don't use the product in a wet ear.
**It may cause damage to the ear.*
- Do not dispose of the batteries by throwing them into a fire.
** They may explode.*

IMPORTANT

- You are advised that changes and modifications not carried out on the equipment by the MINILAND technical service will not be covered by the product guarantee.
- The warranty excludes the user-caused failures listed below:
 - Unauthorized disassembly and modification of the product
 - If the product is dropped during application or transportation
 - An incorrect use of the product.

3.3. RECOMMENDATIONS


- The device is intended for occasional reading and monitoring of body temperature in a domestic environment. Do not use the product for other purposes.
- It is not waterproof. It should only be cleaned with a damp cloth.
- The product should not be exposed to any chemical solvent, direct sunshine or high temperatures.
- Do not drop, step on or subject the product to any vibration or impact.
- Do not use a mobile phone near the thermometer when it is being used.
- For the disposal of batteries and waste, please act according to local law.
- Remove the batteries if the unit is not going to be used for a long period of time.

4. PRODUCT FEATURES

The device is intended for the occasional reading and monitoring of human body temperature in a domestic environment.


4.1. PRODUCT STRUCTURE

4.1.1. PARTS DESCRIPTION

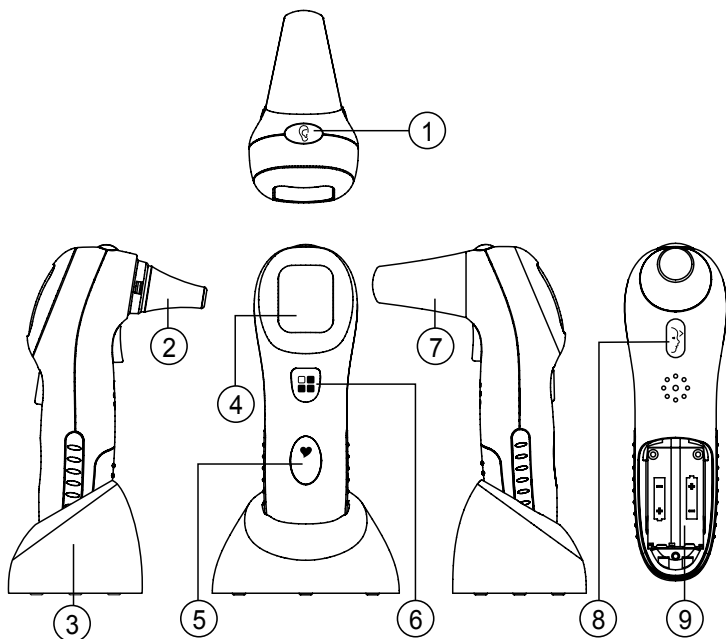
1. Ear key
2. Probe
3. Support
4. LCD
5. Heart rate key
6. Broadcast key
7. Forehead probe
8. Forehead key
9. Battery

4.1.2. LCD DISPLAY

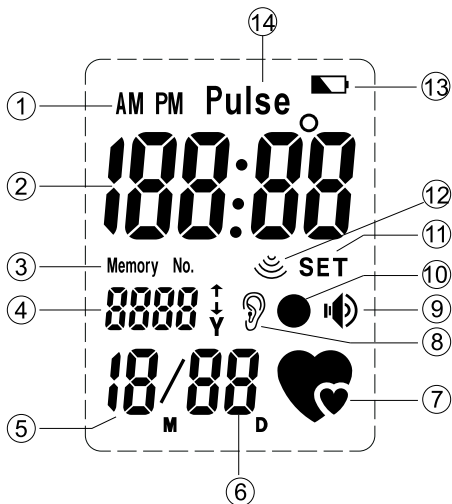1. AM/PM
2. Time, temperature or pulse display
3. Memory number
4. Memory number or year display
5. Month display
6. Date display
7. Heart rate reading symbol
8. Ear temperature reading symbol
9. Speaker
10. Forehead temperature reading symbol
11. Set display
12. Body temperature reading symbol
13. Low battery symbol
14. Pulse

4.2. BASIC FUNCTIONS

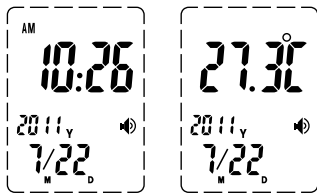1. Environment temperature, current date and time display:

Range of date display falls within Gregorian calendar year 2001~2099, time display mode is 24 hours shift.

LCD can display date, time and environment temperature, LCD displays time and temperature every 5 seconds.

2. Reading and display of temperature:

Readings range from 32.0°C to 42.9°C (89.6°F to 109.3°F).

3. 15 sets of readings:

Maximum of 15 sets of temperature records stored.

Temperature is stored together with date, memory number and mode (Ear or Forehead).

4. Reading display of room temperature :

Reading ranges from 5.0°C to 59.9°C (41.0°F to 139.8°F).

5. Voice prompt: readings, and/or room temperature, and/or current time.

4.3. ADDITIONAL FUNCTION

1. Function: Reading and display of heart rate.

2. Range: 40~140


3. To monitor your heart rate:

If you want to monitor your heart rate, please press [Heart rate] button, the device makes a “DI” sound, the LCD will display as in figure 2 and the device starts to read your heart rate. Meanwhile place your finger on the [Heart rate] key as in figure 1. When the device detects the heart rate, the symbol [♥] will flash to the beat of your heart.

When the reading is complete (about 15 seconds), you will hear a “DI” sound and the LCD will display your heart rate as in figure 3.

Caution:

1. Do not keep the [Heart rate] button pressed down when the device is taking a reading of your heart rate.
2. If the symbol [♥] does not flash, please make sure your finger is placed correctly and you are pressing down hard enough. Excessive pressure will obstruct the local blood circulation and the device will not be able to detect the heart rate.
3. The maximum reading time is 15 seconds.

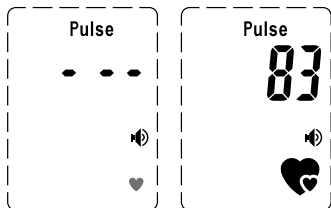5. INSTRUCTIONS FOR USE

5.1. PREPARATIONS

1. Replace any low-voltage battery to ensure full power supply.
2. For a correct reading, the thermometer must be in the same environment 30 minutes before the reading takes place.
3. Any sudden changes in the temperature of the environment are likely to affect an accurate reading.
4. The reading must be taken in a clean ear to assure maximum accuracy.
5. The reading must be taken when the patient is in a state of rest for maximum accuracy.
6. The temperature is likely to rise immediately after exercise or bathing.

5.2. HOW TO TAKE A READING

1. To take an ear temperature:

Pull out the forehead probe from the unit, and insert the temperature probe along the eardrum as deep as possible into the ear. Press [Ⓜ] key to take the ear temperature. Upon completion you will hear two short beep sounds, the LCD will display “ooo” and the

result of the reading, and you will hear the result via the speaker if the TALK function is enabled. When the “ooo” disappears from the screen, the next reading can be taken.

2. To take the forehead temperature:

Put on the forehead probe onto the unit, then locate the probe gently in the center of the forehead, midway between the eyebrow and the hairline. Press [] key once to start the forehead measurement (the measuring time is approximately 1~2 seconds). Upon completion you will hear two short beep sounds, the LCD will display the result of the reading, and you will hear the result via the speaker if the TALK function is enabled.

3. The result will automatically be recorded in the memory. Once 15 recordings have been stored, on taking the 16th recording, the first set of data will be removed and the new recording will take the top position.

4. Three readings can be taken consecutively. However, you should wait for at least 10 minutes before taking a new reading.

5. Display of the results:

A. An accurate value display for result ranges between 32.0°C and 42.9°C (89.6°F and 109.3°F)

B. “HI” is displayed for result above 42.9°C~109.3°F

C. “LO” is displayed for result below 32.0°C~89.6°F

6. Once the reading has finished, the unit will automatically return to time mode after 30 seconds.

7. The display will automatically turn off after 3 minutes after the last key has been pressed (if sleep mode is enabled).

5.3. ADVANCED OPERATION


5.3.1. Clock mode

1. Display time in a sequence of year, month, date, hour, minute with a flickering “:” stands for seconds. The LCD displays the current time and room temperature every 5 seconds.

2. Press [] to show the current time and/or room temperature (except for French and Russian).

3. Press hold [] for 3 seconds to select English, German, French, Spanish, Italian, Russian or the mute function.


You will hear a beep while a reading is being taken if the unit is in talking mode.

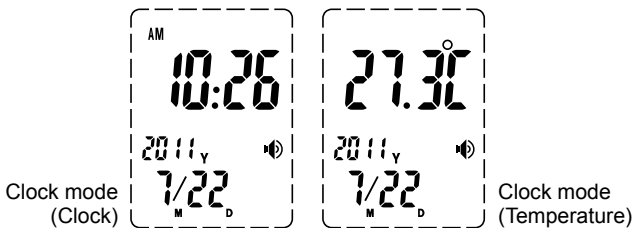4. When the TALK function is enabled, the [] icon is displayed. If it is disabled, the [] icon will be off. The default TALK function is enabled.

5. Press [] and [] key to retrieve the recorded temperature; press the [] and [] key for approximately 3 seconds to enter the set up mode, the time, date and LCD state (on or off) can be reset.

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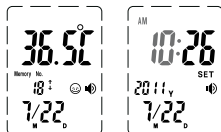5.3.2. SETUP MODE

Press the [] and [] keys for approximately 3 seconds to enter the setup mode.

1. Press [] to select between a setting of HOUR, MINUTE, YEAR, MONTH, DAY and SLEEP enable/ disable and to increase the digit.

2. Press [] to confirm the settings and go to the next setting.

3. Press hold [] to scroll through the digits quickly.

**The unit will automatically return to clock mode 30 seconds after the last key has been pressed.*

Set hour

Set minute

Set year

Set month

Set day

Set sleep
enabled

Set sleep
disabled

5.3.3. MEMORY MODE


1. Press hold [] and [] in clock mode to enter memory mode;

2. “-” will be displayed if no data is stored;

3. Maximum storage capacity is 15 sets of data with date and memory number shown;

4. The unit automatically returns to clock mode 5 seconds after the last key has been pressed.

Nota: Press any key to stop the sound during the announcement.

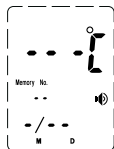Memory no
data

Memory data

5.4. RECOMMENDATIONS

1. When using for the first time, we recommend taking the reading 3 times in the same ear. The maximum reading should be recorded in the case of different results as there may be a slight difference of between (+/- 0.3°C, +/- 0.5°F) if the thermometer is used continuously.

2. Wait only 30 seconds before taking another reading on the same patient again to avoid excessive cooling of the skin.

3. For users in sound health, both ears are almost the same in terms of body temperature. You should select the ear that shows a stable and relatively higher temperature.

4. Please clean the ear before taking a reading.

NOTE: The unit may not function correctly in an environment with electrostatic discharge so it may be necessary to reset the unit.


5.5. TROUBLES AND TROUBLESHOOTING

PROBLEMS	CHECK LIST	SOLUTIONS
No response/ Automatic reset when the insulator is taken out	Battery used up	Insert a new battery
	Battery in wrong polarity	Take out battery, reinsert battery correctly
	Poor battery contact	
Battery symbol on LCD	Low battery life	Insert a new battery
"Lo" on screen (result below 32.0°C or 89.6°F)	Thermometer placed correctly in the eardrum / forehead	Follow user's manual to point at the eardrum / forehead.
"Hi" on screen (result below 42.9°C or 109.3°F)	Please check the operation method	Follow user's manual for correct reading.
Errp" on screen	Hardware problem	Consult with your dealer
"ErrE" on screen	Environment temperature beyond the scope of 10°C~40°C (50.0°F~104.0°F) for Ear 15°C~40°C (59.0°F~104.0°F) for Forehead	Condition the thermometer until it reaches environment temperature 10°C~40°C (50.0°F~104.0°F) for Ear 15°C~40~ (59.0°F~ 104.0°F) for Forehead
"Err" on screen	Operation error, More than 15 seconds without test result in heart rate mode	Please repeat the operation, following the instructions in user manual

6. MAINTENANCE

6.1. REPLACE A BATTERY

Remove the battery compartment flap. Use only new batteries (2 x AAA battery). Replace the battery compartment firmly.

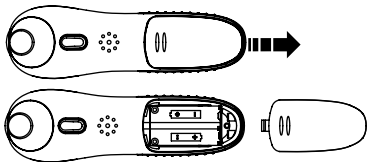6.2. MAINTENANCE, STORAGE AND CALIBRATION

Maintenance:

1. Remove any stains on the casing with a soft dry cloth.
2. Do not wash with water or detergent containing abrasive elements or benzene.

Storage:

1. Clean the thermometer and replace the cover to protect the lens after each use.
2. Do not expose the thermometer to direct sunshine, high temperatures, damp conditions, fire, flame, vibration or impact.
3. Remove the battery for long time of non use.

Calibration:

This thermometer is calibrated at the time of manufacture. If the thermometer is used according to instructions, periodic recalibration is not required. If at any time you question the accuracy of the measurement, please contact the retailer immediately.

Don't attempt to modify or reassemble the thermometer.

6.3. CARE AND CLEANING

1. Keep the probe surface clean, otherwise the reading will not be accurate.
2. Clean the probe with clean paper moistened in water or alcohol, and take the reading after the water on the surface of the probe all evaporates.

Caution:

Make sure the tip of thermometer is clean. A dirty lens may cause an inaccurate reading.

7. TECHNICAL SPECIFICATION

- Product Name: Clinical IR thermometer
- Model: 89068
- Power supply: DC3V (2 x AAA battery)
- Temperature sensing Part: IR sensor
- Range: 32.0°C ~ 42.9°C (89.6°F ~109.3°F)

- Room Temperature: 5.0°C ~ 59.9°C (41.0°F ~139.8°F)
- Accuracy:
+/- 0.2°C (or +/- 0.4°F)
35.5°C ~ 42.0°C (or 95.9°F ~107.6°F)
+/- 0.3°C (+/- 0.5°F)
35.0°C ~ 35.4°C (or 89.6°F ~ 95.8°F)
42.1°C ~ 42.9°C (or 107.7°F ~ 109.3°F)
- Room Temperature: +/- 2°C (4°F)
- Size: 127(L) x 39(W) x 63(D) mm
- Weight : 67g (without batteries)
- Working environment: 10°C~40°C (50.0°F~104.0°F) for Ear
15°C~40°C (59.0°F ~104.0°F) for Forehead
Relative humidity: 20%RH~85%RH
- Storage environment: -10°C~55°C (14°F~131°F)
Relative humidity: 20%RH~85%RH
- Precision/repeatability: +/- 0.3°C (+/- 0.5°F)

**Dispose of used batteries responsibly, please consult your distribution center for details.*

NORMALIZED SYMBOLS

- Attention: Read the operating Instructions for use!

- Applied part of type B

- Disposal in accordance with Directive 2002/96/EC (WEEE)

0197

- Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements. Notified body TÜV Rheinland (CE0197)

- Batch code

- Serial number

- The name and the address of the manufacturer

8.BATTERY & PRODUCT DISPOSAL INFORMATION

- Remember to dispose of the product responsibly.
- Do not mix products bearing the crossed -out bin symbol with your general household waste.
- For the correct collection and treatment of these products take them to the collection points designated by your local authority. Alternatively contact the retailer who sold you the product.

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We declare under our own responsibility that the product:

thermotalk plus · Miniland Baby · 89068 · Miniland, S.A.

to which this declaration refers conforms with the relevant standards or other standardising documents:

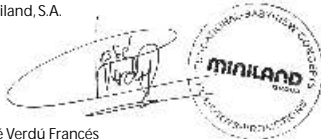
- Directive 93/42/EEC Medical devices
- EN 60601-1:1998+A1:1991+A2:1995 Medical electrical equipment
 - Part 1: General requirements for safety
- EN 60601-1-2: 2007 Medical electrical equipment
 - Part 1-2: General requirements for basic safety and essential performance
 - Collateral standard: Electromagnetic phenomena - Requirements and tests
- EN ISO 13485:2003/AC:2007 Medical devices- Quality management systems
 - requirements for regulatory purposes
- ISO 10993-1 Biological evaluation of medical devices - Part 1: Evaluation and testing
- ISO 10993-5:2009 Biological evaluation of medical devices
 - Part 5: Tests for in vitro cytotoxicity
- ISO 10993-10:2002+Amd 1:2996 Biological evaluation of medical devices
 - Part 10: Tests for irritation and skin sensitization
- ISO 10993-12:2007 Biological evaluation of medical devices
 - Part 12: Sample preparation and reference materials
- EN 12470-5:2003 Clinical thermometers - Part 5: Performance of infra-red ear thermometers (with maximum device)
- Pls delete it, It's no available for our product
- Directive 2006/66/EC, Batteries and accumulators and waste batteries and accumulator (Heavy metals in battery)
- Directive 2002/95/EC RoHS (Restriction of the use of certain hazardous substances in electrical and electronic equipment)
- REACH Regulation (Registration, Evaluation, Authorisation and Restriction of Chemical substances)
- Directive 2002/96/EC Waste electrical and electronic equipment (WEEE), EN 50419:2006 for symbol

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© Miniland, S.A. 2011