miniLAND
BABY

thermoadvanced plus

Castellano	3	Italiano	58
English	14	Česky	69
Português	25	Polski	80
Français	36	Русский	91
Deutsch	47		

www.minilandbaby.com

ÍNDICE

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INDEX

1. Introduction
2. Contents
3. Safety instructions
4. Product features
5. Instructions for use
6. Maintenance
7. Technical specification
8. Battery & product disposal information

miniLAMP

English

1. INTRODUCTION

Congratulations on the purchase of this IR clinical speaking thermometer

Before using this device, please read all the instructions in this manual carefully. Please keep in a safe place as it contains important information.

The features described in this user's guide are subject to modifications without prior notice.

2. CONTENTS

- 1 IR no contact speaking thermometer
- 1 storage and travel bag
- 2 AAA Batteries
- 1 Instruction manual and guarantee

If any of the above is missing, please contact your retailer.

3. SAFETY INSTRUCTIONS

- 1.This device is not a toy; do not allow children to play with it.
- 2.Keep the thermometer out of the reach of children.
- 3.Read the instructions carefully before use and save them for future reference.
- 4.Use of this thermometer is not intended as a substitution for consultation with your physician. Reading results are for reference only. Contact your physician if you have or suspect any health concerns.
- 5.Make sure the batteries are inserted with the correct polarity.
- 6.Consumed batteries should be removed.
- 7.Don't throw battery into a fire, the battery is likely to explode.
- 8.Never attempt to open any part or dismantle the controls. Responsibility for failure cannot be accepted if these instructions have been ignored.
- 9.Don't take apart, repair or change any parts of the unit at any time.
- 10.Don't use the thermometer if the main body is damaged (for example, the infrared body is broken). The continuous use of a damaged unit may cause injury, erroneous results or serious danger.
- 11.Do not shake or tap the thermometer. Do not drop it.
- 12.The infrared sensor must be kept clean, dry and undamaged all the times to ensure accurate readings.

13. Do not operate the system in the vicinity of strong transmitter such as a microwave oven and shortwave appliances.

14. This product needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided, and this unit can be affected by portable and mobile RF communications equipment.

15. To ensure eMyBaby App readings are accurate, keep at a minimum distance of about 30 cm from other wireless devices.

16. Don't store the thermometer in extreme temperatures below -20°C/-4°F or over 55°C/131°F) or extreme humidity (below 10%RH or over 90%RH). This could cause inaccurate readings.

17. This device is intended for occasional reading and monitoring of the human body or objects temperature in a domestic environment. Do not use the product for other purposes.


You are advised that changes and modifications not carried out on the equipment by the MINILAND technical service are not covered by the product guarantee.

4. PRODUCT FEATURES

4.1. PARTS DESCRIPTION


4.1.1. THERMOMETER

1. Forehead temperature measurement button
2. Infrared sensor
3. Spotlight
4. Battery cover
5. Set/Memory button
6. Loudspeaking
7. Talking button
8. Object temperature measurement button
9. LCD

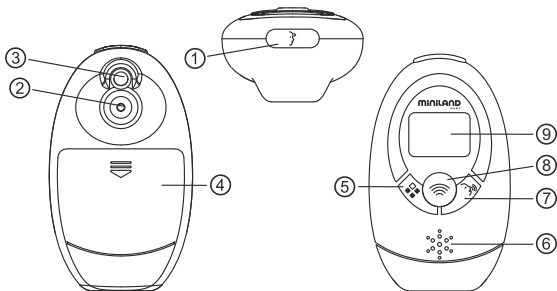4.1.2. LCD DISPLAY

- A. Temperature display
- B. PM (time format)
- C. Memory sign
- D. Low battery sign
- E. Room temperature measurement sign
- F. Object temperature measurement sign
- G. Forehead temperature measurement sign

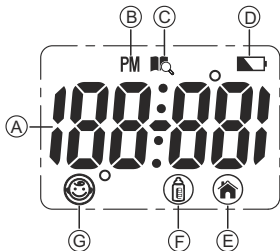5. INSTRUCTIONS FOR USE

5.1. BEFORE USE

5.1.1. RECOMMENDATIONS

1. To achieve an accurate reading, it is very important to check the infrared sensor is clean before using.
2. If the temperature where the thermometer has been stored differs greatly from where the temperature is going to be taken, please wait until the thermometer temperature stabilizes to the room temperature, about 30 minutes before using.
3. Don't touch or blow the infrared sensor. A dirty infrared sensor may cause an inaccurate reading.
4. Remember that in waiting mode, the LCD alternatively displays time, date and room temperature.

5.1.2. INSTALLATION/REMOVAL OF BATTERIES

When the low battery indicator (D) appears on the LCD, or the thermometer does not function at all, you should replace the old batteries with new ones as soon as possible.

This product uses two AAA 1.5V batteries, use the following procedure to replace or insert batteries:

1. Open the battery cover following the direction of the arrow.
2. Remove the old batteries and then insert the new batteries into the battery compartment. Make sure you install the new batteries properly; be sure to match the positive and negative ends of the battery to the marking in the batteries compartment of the thermometer.
3. Replace the battery compartment cover and press to close.

5.2. FOREHEAD TEMPERATURE MEASUREMENT

- Make sure the probe is clean before measuring.
- Make sure that the patient's forehead is clean before taking a reading to assure that it is accurate. If it is necessary to wipe the patient's forehead, please wait for 5-10 minutes before taking his temperature.
- External factors may affect the patient's forehead temperature, particularly when an individual has been exposed to very hot or very cold temperatures or has recently been swimming or bathing. If necessary, wait 30 minutes before taking the patient's temperature.
- A forehead and temple temperature differs from an internal temperature, which

is taken orally or rectally. Vasoconstriction, an effect which constricts the blood vessels and cools the skin, can occur during the early stages of a fever. In this case, the temperature taken may be unusually low. If the reading does not match the patient's perception of his own body temperature, repeat the reading every 15 minutes. As a reference point, you can also measure the internal body temperature using a conventional oral or rectal thermometer.

·The baby's skin reacts very quickly to the ambient temperature. Therefore, do not take the baby's temperature with the thermometer during or after breastfeeding, because her skin temperature may be lower than the internal body temperature.

1. Place the thermometer in your hand with your finger in the forehead button (1).

2. Hold the thermometer above the center of the forehead, at a distance approximately 3-5 cm, midway between eyebrow and hairline.



3. Press and hold the forehead button (1), the spotlight is ON and you can locate the correct position for a forehead reading. Hold the forehead button as long as you want the spotlight to stay on. If the position is correct, release the forehead button to start the reading and then the thermometer beeps once.

4. Keep the thermometer flat until you hear two beeps and the temperature is given. Then the backlight will be lit up according to temperature taken.

5. The reading is complete and the thermometer may be removed from the head.

6. Read the temperature on the display

7. You can press the forehead button (1) to start a new reading or press Set/Mem button (5) to return to waiting mode.

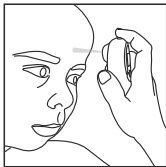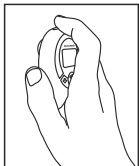NOTE

The spotlight will stay on whilst the reading button is being pressed. The spotlight turns off after you release the button. The spotlight will automatically turn off if you press the forehead button (1) for more than 1 minute; after you release the button, the device reverts back to waiting mode.

Normal body temperature is within a set range. A person's normal temperature range tends to decrease with age. The following table shows a normal temperature ranges for age:

0-2 years	36.4°C-38.0°C	97.5°F-100.4°F
3-10 years	36.1°C-37.8°C	97.0°F-100.0°F
11-65 years	35.9°C-37.6°C	96.6°F-99.7°F
> 65 years	35.8°C-37.5°C	96.4°F-99.5°F

The range varies from person to person and can be influenced by many factors such as the time of day, the level of activity, medication, emotion and so on. So we recommend that you practice with the thermometer on yourself and family members when you are healthy. This way you know how the thermometer works and can feel more confident of the readings you take when a family member is ill.

NOTE: If the measured temperature is equal or greater than 38°C, the display will show the fever indicator and will be backlit in red.

5.2.1. TEMPERATURE BACKLIGHT INDICATION

The display backlight will light up after a human body temperature reading has finished, as follows:

Blue backlight: below 37.9°C (100.3°F)

Red backlight: 38.0°C-42.9°C (100.4°F-109.3°F)

5.3. OBJECT TEMPERATURE READING

- Make sure the probe is clean before taking the temperature.
- Don't use this mode to take someone's body temperature.
- When reading the temperature of an object, the current surface temperature is displayed. It can be different from the internal temperature, especially if the surface is exposed to direct sunlight or a draught.

1. Place the thermometer in your hand with your finger in the object button (8).

2. Hold the thermometer so that the probe position is at a distance within 3-5 cm from the liquid and/or surface you want to take the temperature of.

3. Press the object button (8) to start the reading, and then the thermometer beeps once.


4. Remove the thermometer from the object. The LCD displays the temperature.

5. You can press the object button (8) to start a new reading or press Set/Mem button (5) to return to waiting mode.

5.4. CHECK THE MEMORY

There are a total of 9 memory data for recording forehead temperature readings. The current reading is always stored in the last storage space. When all storage spaces have been occupied, the oldest reading is always deleted from the memory.

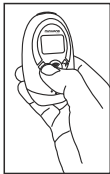1. Press Set/Mem button (5) to check the memory number, temperature and date.
2. Press Set/Mem button (5) to check the next temperature of the memory.
3. Press and hold Set/Mem button (5) to check all the memory data quickly.
4. Press object (8) and forehead (1) measuring button at the same time, the device will revert back to waiting mode.

5.5 DATE AND TIME SETTING

1. Press and hold Set/Mem button (5) about 2 seconds to enter into setting mode under the waiting mode.

2. Press again Set/Mem button (5) to select between a setting of hour mode (12hr or 24hr, "PM" will display on LCD when you select 12hr mode), hour, minute, year (2001-2099), month, day, year and sleep mode (ON or OFF).

3. Press talk button (7) to confirm the settings and go to the next setting.

When the sleep mode is activated, the thermometer will turn off automatically after approximately 3 min; when the sleep mode is deactivated, the thermometer remains power on. To save power, we recommend to set sleep mode on.


5.6. VOICE PROMPT

In waiting mode, press the talk button (7) to report the temperature, current time or room temperature.

After each forehead reading, the thermometer reports the temperature obtained in the language chosen.

This thermometer communicates the temperature with a voice message, 6 languages: Spanish, English, German, French, Italian and Russian.

· To choose the desired language or mute option, in waiting mode press and hold talk button (7) about 2 seconds to select the talking language which you want, it will change to another language or mute function automatically.

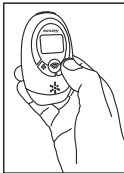5.7. SWITCHING BETWEEN FAHRENHEIT OR CELSIUS When the LCD displays a reading or memory temperature, press and hold talk button (7) about 3 seconds this will change the mode to either °F or °C.

5.8. BABY TEMPERATURE CONTROL

This thermometer has a special communications module that lets you connect to your smartphone or tablet, automatically sending the results of each reading, to keep track of temperatures, generate graphs and send them by e-mail.

5.8.1. RECOMMENDATIONS

1. Make sure that your smartphone or is running on the latest operating system. For information about software updates, see the instructions for your device.
2. Make sure that your eMyBaby software (app) is updated at regular intervals. For information about app updates, see the instructions for your device.
3. For the communications module to function correctly, the mobile device has to be sufficiently charged. Users should have a power supply unit or a charging device to hand.
4. Make sure that there is always enough memory available to back up the reading data on your smartphone or tablet.
5. If you receive a call or a text message during a reading, the reading could be cancelled or it can provide an inaccurate information. We recommend switching your device to flight mode.

CAUTION:

This machine should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, this machine should be observed to verify normal operation in the configuration in which will be used.

5.8.2. eMyBaby App

Connect with smartphone or tablet for recording and controlling baby temperatures.

To operate your new infrared thermometer with smartphone or tablet, you need to install the free eMyBaby App.

Once the software has been successfully installed, the new eMyBaby icon  appears on your device.

Log in the eMybaby application, select the option "thermometer" in the main menu and follow the instructions listed there.

NOTE

The connection to a mobile device may be not always be successful in noisy surroundings.

Connection with Apple mobile devices: iPhone, iPad & iPod touch (iOS 4.3 or later)

5.9. CALIBRATION

The thermometer is initially calibrated at the time of manufacture. If the thermometer is used according to the instructions, periodic recalibration is not required. If any time you question the accuracy of the reading, please contact the retailer immediately.

Don't attempt to modify or reassemble the thermometer.

5.10. TROUBLESHOOTING

PROBLEMS	CHECKLIST	SOLUTIONS
No response/Automatic reset when pull out insulator	Battery used up	Change new battery
	Battery in wrong polarity	Take out battery, reinsert battery correctly.
	Poor battery contact	
"Lo" on screen	The measured temperature is lower than 32.0°C/89.6°F. Please check the operation method.	Follow the steps in User's Manual for proper reading.
"Hi" on screen	The measured temperature is higher than 42.9°C/109.3°F. Please check the operation method.	Follow the steps in User's Manual for proper reading.
"ErrP" or "ErrH" on screen	Hardware problem	Contact the manufacturer or distributor
"ErrE" on screen	Environment temperature beyond the scope of 10°C~40°C (50.0°F~104.0°F) for Objects 15°C~40°C (59.0°F~104.0°F) for Forehead	Condition the thermometer until it reaches environment temperature 10°C~40°C (50.0°F~104.0°F) for Objects 15°C~40~ (59.0°F~104.0°F) for Forehead
The reading is not accurate or if there is any doubt with the result	Please check if the infrared sensor is clean or not	Clean the infrared sensor with cotton swab according to user manual.
	Please check if the reading mode is correct or not	Ensure you have read the manual and know how to use the thermometer properly.
	Please check if you have let the thermometer and patient stabilize in the room for at least 30 minutes	Please keep the thermometer and patient in the room where the temperature is to be taken at least 30 minutes before using
	Are you using the thermometer indoor	Please take the reading indoor
Battery symbol on LCD	Low battery life	Replace with new battery
Other malfunctions	Contact the MINILAND technical service	

6. MAINTENANCE

- Please place away from a draught to assure an accurate reading.
- Prevent the product getting wet, if it does get wet please dry thoroughly. Please note that corrosive liquids may cause damage.
- Please use and preserve at a normal temperature, the product life will be shortened if temperatures are excessively high or low.
- If your thermometer will not be used on a regular basis, remove the battery to prevent possible damage due to chemical leakage. If battery does leak, remove carefully. Do not allow bare skin to come in contact with the leaking fluid.

6.1. CARE AND CLEANING

- To ensure accurate measurements, keep the probe tip clean and free of scratches. Finger prints or dirt will affect the accuracy of the thermometer.
 - In order to get an accurate reading, please clean the probe tip.
 - To clean the sensor, gently wipe its surface with a cotton swab slightly moistened with alcohol and immediately wipe dry with a clean cotton swab.
 - Use a soft dry cloth to clean the body of the thermometer. Never clean the thermometer with abrasive cleanser, thinner, benzene or submerge the thermometer into water or other liquids.
- After cleaning, allow at least 20 minutes drying time before taking temperatures. Place the thermometer in the protective case or table stand holder. Store it in a dry, clean place at room temperature

7. TECHNICAL SPECIFICATION

Power supply	DC 3V (2 x 1.5 AAA batteries)	
Body reading range	32.0°C-42.9°C	89.6°F-109.3°F
Forehead reading accuracy	In range of 35.5°C-42.0°C 95.9°F-107.6°F	± 0,2 °C ± 0,4 °F
	For other range	± 0,3 °C ± 0,5 °F
Object measure range	20,0 °C-60,0 °C	68,0 °F-140,0 °F
Object measure accuracy	± 1,0 °C	± 2,0 °F
Room temperature range	5,0 °C-59,9 °C	41,0 °F-139,8 °F
Room temperature accuracy	± 1,0 °C	± 2,0 °F
Dimension	89 (L) x 53 (A) x 35 (F) mm	
Weight	56 g (without battery)	
Forehead operation condition	With relative humidity of 20%-85%	
	15,0 °C-40,0 °C	59,0 °F-104,0 °F

Object operation condition	With relative humidity 20%-85%	
	10,0 °C-40,0 °C	50,0 °F-104,0 °F
Storage environment	With relative humidity 20%-90%	
	-20,0 °C-55 °C	-4,0 °F-131,0 °F

7.1. NORMALIZED SYMBOLS

Type BF Applied Part.

Disposal in accordance with Directive 2002/96/EC (WEEE).

Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements.

0197

Notified body TÜV Rheinland (0197).

The name and the address of the manufacturer.

The name and the address of the Authorized EC-representative in Europe..

Refer to instruction manual.

8. BATTERY & PRODUCT DISPOSAL INFORMATION

Remember to dispose of the product responsibly.

Do not mix products bearing the crossed out bin symbol with your general household waste.

For the correct collection and treatment of these products take them to the collection points designated by your local authority. Alternatively contact the retailer who sold you the product.

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We declare under our own responsibility that the product:


thermoadvanced plus · Miniland Baby · 89083 · Miniland, S.A.

**to which this declaration refers conforms with the
relevant standards or other standardising documents:**

- Directive 93/42/EEC Medical devices and amended by Directive 2007/47/EC (CE0197). (EC certification No: 170 198 16001)
- EN 60601-1:2006 Medical electrical equipment - Part 1: General requirements for safety
- EN 60601-1-2:2007 Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral standard: Electromagnetic phenomena - Requirements and tests
- EN ISO 13485:2012 Medical devices- Quality management systems - requirements for regulatory purposes
- ISO 10993-5:2009 Biological evaluation of medical devices - Part 5: Tests for in vitro cytotoxicity (Report No: SBRTC-2010-126)
- ISO 10993-10:2002+Amd 1:2996 Biological evaluation of medical devices - Part 10: Tests for irritation and delayed-type hypersensitivity AMENDMENT 1
- EN 12470-5:2003 Clinical thermometers - Part 5: Performance of infra-red ear thermometers (with maximum device)
- ISO 14971 :2012: Medical devices -Application of risk management to medical devices
- Directive 2006/66/EC, Batteries and accumulators and waste batteries and accumulator (Heavy metals in battery)
- Directive 2002/95/EC RoHS (Restriction of the use of certain hazardous substances in electrical and electronic equipment)
- REACH Regulation (Registration, Evaluation, Authorisation and Restriction of Chemical substances)
- Directive 2002/96/EC Waste electrical and electronic equipment (WEEE), EN 50419:2006 for symbol

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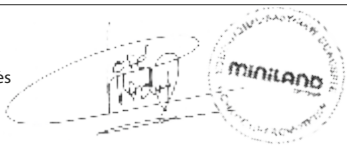Made in China. Keep this information.

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Selective pick-up. Protect the environment.

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© **Miniland, S.A. 2013**
Miniland S.A. P.Ind. La Marjal I C/ La Patronal, 8-10
Shenzhen Dongdixin Technology Co., Ltd
Shanghai International Holding Corp. GmbH (Europe)